

# PlancheHisto Manual

28/06/2022

sebastien.schaub@imev-mer.fr

## I. Goal

The program helps to adapt white balance between images to get homogeneous gallery. It uses the ActionBar from J.Mutterer

## II. Installation


Program, tutorial are available here : <https://github.com/SebastienSchaub/PlancheHisto>

A video tutorial is available here : <https://twitter.com/i/status/1409864168749027329>

- For simplicity the program and icon images are incorporated in the directory of ImageJ.
- Copy the folder "macros" and "plugins" in the Fiji/ImageJ folder (no file will be overwritten except action\_bar.jar if you already have).
- Open ImageJ
- Plugins > Macros > Run... and select [Home Directory]\macros\MICA\PlancheHisto.ijm.

## III. How to Use

Here are the functions linked to the program :

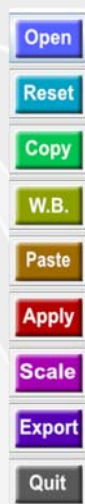

**[1] Open :** To open an image. It convert the image in 3 channels (R,G,B).

**[2] Reset :** To reset the white balance per channel from 0 to intensity max.

**[3] Copy :** copy the ratio R,G,B range to apply to other image

**[4] W.B. :** Apply automatic white balance on the current ROI (if no ROI, apply it on the whole image)

**[5] Paste :** Apply white balance on the image to make that the ROI has the same ration R,G,B as in the copied ROI.

**[6] Apply :** Apply the same transformation to the current image.

**[7] Scale :** Show or hide the scale bar

**[8] Export :** Export the image in RGB color coded image

**[0] Parameter :** give access to parameter of the image and let scale the image (if not done):



**Institut de la Mer de Villefranche/Mer**


181, Chemin du Lazaret,

06234 Villefranche-sur-Mer Cedex, France


Sebastien.Schaub@imev-mer.fr



## Plateforme d'Imagerie par Microscopie (PIM)

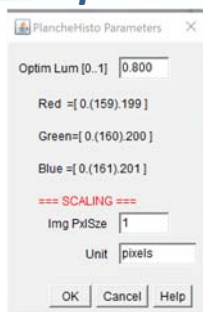

- Optim Lum [0..1]: is the brightness which is reached with white balance. The highest, the brightest
- Red : gives intensities : [Minimum, Modal, Maximum]
- Green : idem
- Blue : idem
- ImgPxSize : the pixel size measured in [Unit]
- Unit : pixel size scale [pixel/microns/ etc....]
- [Help] : provides information to calibrate images from PIM Microscopy platform (feel free to change it based on your needs).

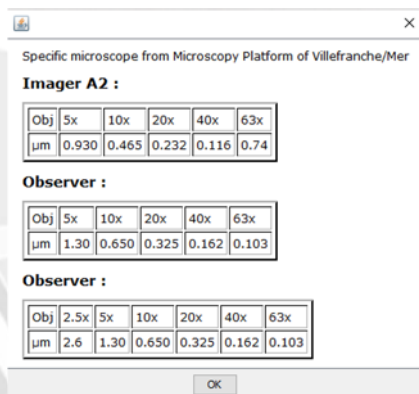

**Institut de la Mer de Villefranche/Mer**

181, Chemin du Lazaret,

06234 Villefranche-sur-Mer Cedex, France



Sebastien.Schaub@imev-mer.fr